INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

C12N 5/08, A61K 121/00	(51) International Patent Classification ⁶ :		(11) International Publication Number: WO 95/20649
(22) International Filing Date: 6 January 1995 (06.01.95) (30) Priority Data: 08/188,262 27 January 1994 (27.01.94) 08/214,400 16 March 1994 (16.03.94) 08/300,982 6 September 1994 (06.09.94) (71) Applicant: CELLCOR, INC. [US/US]; 200 Wells Avenue, Newton, MA 02159 (US). (72) Inventors: BABBITT, Bruce, P.; 142 Massapoag Avenue, N. Easton, MA 02356 (US). ZHANG, Zhengyi, J.; 31 Pheasant Road, Needham, MA 02192 (US). OSBAND, Michael, E.; 101 Blake Road, Brookline, MA 02146 (US). (74) Agent: CLARK, Paul, T.; Fish & Richardson P.C., 225 Franklin	C12N 5/08, A61K 121/00	A1	(43) International Publication Date: 3 August 1995 (03.08.95)
08/188,262 27 January 1994 (27.01.94) US 08/214,400 16 March 1994 (16.03.94) US 08/300,982 6 September 1994 (06.09.94) US (71) Applicant: CELLCOR, INC. [US/US]; 200 Wells Avenue, Newton, MA 02159 (US). (72) Inventors: BABBITT, Bruce, P.; 142 Massapoag Avenue, N. Easton, MA 02356 (US). ZHANG, Zhengyi, J.; 31 Pheasant Road, Needham, MA 02192 (US). OSBAND, Michael, E.; 101 Blake Road, Brookline, MA 02146 (US). (74) Agent: CLARK, Paul, T.; Fish & Richardson P.C., 225 Franklin	(22)		BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL,
	08/188,262 27 January 1994 (27.01.94) 08/214,400 16 March 1994 (16.03.94) 08/300,982 6 September 1994 (06.09.94) (71) Applicant: CELLCOR, INC. [US/US]; 200 Wells Newton, MA 02159 (US). (72) Inventors: BABBITT, Bruce, P.; 142 Massapoag Av Easton, MA 02356 (US). ZHANG, Zhengyi, J.; 31 Road, Needham, MA 02192 (US). OSBAND, Mid 101 Blake Road, Brookline, MA 02146 (US). (74) Agent: CLARK, Paul, T.; Fish & Richardson P.C., 225	Avenu Avenu venue, Pheasa	With international search report. We have the contract of the

(57) Abstract

Disclosed is a process of activating patient-derived mononuclear cells by exposing the cells *in vitro* to substances to generate immunoreactive cells. The *ex vivo* activated cells are then reinfused into the patient to enhance the immune system to treat various forms of cancer, infectious diseases, autoimmune diseases or immune deficiency diseases.

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GB	United Kingdom	MR	Mauritania
ΑÜ	Australia	GE	Georgia	MW	Malawi
BB	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Greece	NL	Netherlands
BF	Burkina Faso	HU	Hungary	NO	Norway
BG	Bulgaria	IE	Ireland	NZ	New Zealand
BJ	Benin	IT	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY	Belarus	KE	Kenya	RO	Romania
CA	Canada	KG	Kyrgystan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic	SD	Sudan
CG	Congo		of Korea	SE	Sweden
СН	Switzerland	KR	Republic of Korea	SI	Slovenia
CI	Côte d'Ivoire	KZ	Kazakhstan	SK	Slovakia
CM	Cameroon	LI	Liechtenstein	SN	Senegal
CN	China	LK	Sri Lanka	TD	Chad
CS	Czechoslovakia	LU	Luxembourg	TG	Togo
CZ	Czech Republic	LV	Latvia	TJ	Tajikistan
DE	Germany	MC	Monaco	TT	Trinidad and Tobago
DK	Denmark	MD	Republic of Moldova	UA	Ukraine
ES	Spain	MG	Madagascar	US	United States of America
FI	Finland	ML	Mali	UZ	Uzbekistan
FR	France	MN	Mongolia	VN	Viet Nam
GA	Gabon		· ·		

- 1 -

EX VIVO ACTIVATION OF IMMUNE CELLS

Background of the Invention

5 This invention relates to immunotherapy.

Adoptive immunotherapy as an approach to treating cancer has been evolving rapidly over the past two decades. Beginning in the 1980's, IL-2 and lymphokine-activated killer (LAK) cell-based immunotherapies were shown to induce tumor responses in patients with melanoma and renal cell carcinoma (Rosenberg et al., 1985, N. Eng. J. Med. 313:1485-1492; Rosenberg et al., 1986, Important Advances in Oncology, DeVita et al., eds, Philadelphia, Lipincott, p. 55; Schoof et al., 1988, Cancer Res.

- 15 48:5007-5010; West et al., 1987, N. Engl. J. Med. 316:898-905). However, these therapies require systemic administration of high dose IL-2 which usually is accompanied by severe toxicity (Siegel et al., 1991, J. Clin. Oncol. 9:694-704; Margolin et al., 1989, J. Clin.
- 20 Oncol. 7:486-498). More recently, tumor infiltrating lymphocytes (TIL) isolated from surgically removed tumors were utilized for adoptive immunotherapy following ex vivo activation and expansion (Rosenberg et al., 1986, Science 223:1318-1321; Kradin et al., 1989, Lancet 1:577-
- 25 580; Speiss et al., 1987, J. Natl. Cancer Instit.
 79:1067-1075; Rosenberg et al., 1988, N. Engl. J. Med.
 319:1676-1688). Although treatment of patients with TIL
 has been shown to achieve somewhat higher tumor response
 rates relative to LAK therapy for melanoma (Rosenberg et
- 30 al., 1986, Science 223:1318-1321; Kradin et al., 1989,
 Lancet 1:577-580; Spiess et al., 1987, J. Natl. Cancer
 Instit. 79:1067-1075; Rosenberg et al., 1988, N. Engl. J.
 Med. 319:1676-1688), this type of therapy has not shown
 consistent outcomes in clinical trials of other tumors.
- 35 Furthermore, the limited availability of sufficient

- 2 -

numbers of lymphocytes from autologous tumor tissue and other problems associated with long term, high volume cell culture have also restricted the application of TIL therapy.

5

Summary of the Invention

The invention provides a novel, safe, and costeffective method of nonspecifically enhancing a cellmediated immune response to specific antigens or foreign

10 substances in the body, including cancer cells. The
process involves removing a patient's mononuclear cells
and exposing the cells in vitro to substances which
enhance the immune function of the cells. The ex vivo
activated (EVA) cells are then reinfused into the patient

15 to enhance the patient's immune responses and to treat
various forms of cancer, infectious diseases, autoimmune
diseases, or immune deficiency diseases.

The invention features a process of producing a population of immunoreactive cells by (a) contacting a sample of mononuclear cells derived from a patient, e.g., peripheral blood mononuclear cells (PBMC), with OKT3 at or below 37°C to produce an OKT3-derived culture supernatant (T3CS); (b) removing the T3CS from the sample of patient-derived mononuclear cells; (c) determining the concentration of OKT3 in the T3CS, and if required, supplementing the T3CS with additional OKT3 to achieve a concentration of at least 0.1 ng/ml; (d) providing a second sample of mononuclear cells derived from the patient; and, (e) contacting the second sample of cells with the previously-generated T3CS for a period of time sufficient to yield a population of immunoreactive cells.

A sample of patient-derived mononuclear cells may contain T cells, B cells, monocytes and macrophages as well as other immune cells such as polymorphonuclear leukocytes, neutrophils, eosinophils, natural killer

- 3 -

cells, and stem cells. Mononuclear cells may be derived from the peripheral blood of the patient, or other sites, e.g., a tumor or tumor-draining lymph node.

T3CS is a conditioned medium containing a mixture

of autologous cytokines together with OKT3. The
autologous cytokines mixture preferably promotes the
growth and differentiation of Th1-type T cells, rather
than Th2-type T cells. The OKT3 which is preferably in
solution phase catalyzes the polyclonal activation of

T cells, while the cytokines act synergistically as costimulants to optimize the overall degree of activation.
The presence of both OKT3 and cytokines prevents the
generation of T cells that are anergic or apoptotic and
overcomes signal transduction defects in mononuclear

cells derived from patients with cancer or chronic
infectious diseases.

By the term "immunoreactive cells" is meant polyclonal T cells that exist in a primed state of activation. Primed cells are multifunctional, i.e., they possess an enhanced capacity to proliferate and produce cytokines upon further stimulation. The primed state of activation of the immunoreactive cells induced by culture in the OKT3-autologous cytokine mixture can be identified by measuring the stable biochemical changes, e.g., expression of growth, differentiation, and activation markers, which occur both on the cell surface and intracellularly. Immunoreactive cells of the invention have enhanced immunologic effector function, e.g., helper activity (CD4+ T cells) or cytotoxicity (CD8+ T cells), compared to unprocessed patient-derived mononuclear cells.

Immunoreactive cells have a low spontaneous level of immune function following processing, but are highly sensitized to respond to low doses of second signals upon further culture, or in vivo. The immunoreactive cells of

the invention therefore require further exposure to an immune stimulant, such as an antigen; target cell, e.g., a tumor cell or virus-infected cell; an inflammatory molecule; an adhesion molecule; an immune cell, e.g., an accessory cell; a cytokine; or any combination thereof, to achieve full immunologic effector function. The immunoreactive cells of the invention are multifunctional, polyclonally-activated T cells which have been generated independent of disease-specific antigens utilizing a mixture of nonspecific lymphocyte activators, i.e., autologous cytokines, and a mouse monoclonal antibody, i.e, OKT3, as synergistic stimulants.

Suppressor cells in a population of patient
derived mononuclear cells may be inactivated by contacting the second sample with a suppressor cell inhibitory compound, e.g., cimetidine, indomethacin, cyclophosphamide, ranitidine, pepsid, or any combination thereof. Other histamine type-2 receptor blockers may also be used alone or in combination with the compounds listed above. Cimetidine and indomethacin are preferably used together at concentrations of 5 × 10⁻⁵ M (± 2-fold) and 0.8 × 10⁻⁸ M (± 2-fold), respectively.

The concentration of OKT3 used in the process of
the invention is an amount that promotes activation of
the patient's cells, but leaves minimal surface-bound
OKT3 on the activated cell product. Minimizing the
amount of surface-bound OKT3 on the immunoreactive cells
in turn minimizes human anti-mouse antigen (HAMA) immune
responses and rapid clearance of the immunoreactive cells
from the circulation of a patient undergoing therapy with
the immunoreactive cells of the invention. The
concentration of OKT3 is preferably greater than 0.1
ng/ml but less than 25 ng/ml, more preferably 1-25 ng/ml,
and most preferably 10-15 ng/ml. In addition to OKT3,

- 5 **-**

any compound that binds to the T cell receptor or the T cell receptor-associated CD3 molecule on the cell surface may be used to stimulate the first sample of patient-derived mononuclear cells to produce T3CS.

out for a period of time sufficient to produce a mixture of nonspecific lymphocyte activators capable of promoting the OKT3-catalyzed activation and differentiation of the second sample of patient-derived mononuclear cells into a population of immunoreactive cells. The culture period may range from 1 to 7 days and is preferably 3 days. The length of culture of the first sample may be adjusted, e.g., prolonged, to achieve the desired concentration of a nonspecific lymphocyte activator, e.g., a cytokine, e.g., tumor necrosis factor-alpha (TNFα) or interleukin-2 (IL-2), in the T3CS.

Culture of the second sample of patient-derived cells with T3CS is carried out for a period of time sufficient to produce a population of immunoreactive cells. As discussed above, the immunoreactive cells are in a primed state of activation, i.e., the cells are no longer in a resting state but require an additional stimulus, e.g., exposure to an antigen or other immune stimuli, to achieve a fully activated state characterized by enhanced immune function compared to unprocessed cells. Full activation of immunoreactive cells may be measured by expression of new cell surface markers, e.g., CD25, secretion of lymphokines, e.g., IFNγ, GM-CSF, or TNFα, cellular proliferation, or cellular differentiation into effector cells, e.g., cytolytic T cells or helper T cells.

Culture of patient-derived mononuclear cells with T3CS may be carried out for a period of 1 to 30 days, preferably 5 days. In the absence of antigen, the length of the culture period may be as short as 1-3 days; in the

presence of antigen, the cells may be co-cultured for a period of 1 to 30 days. The length of culture may be adjusted, e.g., prolonged, to achieve the desired level of activation of the immunoreactive cells.

The invention also features a process of 5 producing a population of immunoreactive cells by (a) providing a first sample of mononuclear cells derived from a patient; (b) determining the concentration of Fcreceptor positive accessory cells in the first sample, 10 and if required, supplementing the first sample with a second sample of mononuclear cells derived from the same patient to achieve a concentration of 0.1-50% Fc-receptor positive accessory cells in the first sample; (c) contacting the first sample with OKT3 at or below 37°C to 15 produce a T3CS; (d) removing the T3CS from the first sample; (e) determining the concentration of OKT3 in the T3CS, and if required, supplementing the T3CS with OKT3 to achieve a concentration of at least 0.1 ng/ml; (f) providing a third sample of mononuclear cells derived 20 from the same patient; (g) contacting the third sample with T3CS for a period of time sufficient to activate the third sample in vitro to yield a population of immunoreactive cells. The Fc-receptor positive cells are preferably monocytes, but may be granulocytes or 25 dendritic cells. The concentration of patient-derived monocytes is preferably 0.1-50%, more preferably 1-30%, more preferably 5-15%, and most preferably 10% of the cells in the sample. The second sample of patientderived mononuclear cells may be enriched for monocytes 30 using cell fractionation techniques known in the art, e.g. panning or FACS, prior to augmenting the concentration of monocytes in the first sample.

In another aspect, the process is carried out by:
(a) contacting a first sample of mononuclear cells
35 derived from a patient with OKT3 at or below 37°C to

produce a T3CS; (b) removing the T3CS from the first sample; (c) determining the concentration of tumor necrosis factor-alpha (TNFα) in the T3CS, and if required, supplementing the T3CS with $TNF\alpha$ to achieve a 5 concentration of at least 5 pg/ml; (d) providing a second sample of mononuclear cells derived from the same patient; (e) inactivating suppressor cells in the second sample; and, (f) contacting the second sample with T3CS and for a period of time sufficient to activate the 10 second sample in vitro to yield a population of immunoreactive cells. The concentration of TNF α is preferably in the range of 100 pg/ml to 100 ng/ml, more preferably in the range of 100 pg/ml to 3000 pg/ml, and most preferably in the range of 500 pg/ml to 1000 pg/ml. 15 The concentration of a cytokine, e.g., $TNF\alpha$, may be augmented by (1) altering the length of the T3CSgeneration step, (2) supplementing T3CS with purified non-recombinant cytokine, (3) supplementing T3CS with a recombinant cytokine, or (4) increasing the concentration 20 of a cytokine-producing patient-derived mononuclear cell. The T3CS generation step (step (a)) of the process may be carried out for 1-7 days and is preferably carried out for 3 days.

In yet another aspect, the process requires after step (b), determining the concentration of both TNF α and OKT3, and if required, supplementing the T3CS with TNF α to achieve a concentration of at least 5 pg/ml TNF α and with OKT3 to achieve a concentration of at least 0.1 ng/ml OKT3.

The levels of interferon-gamma (IFN-γ), interleukin-4 (IL-4), and interleukin-10 (IL-10) in T3CS may be adjusted to achieve optimal generation of a particular type of immunoreactive cell, e.g., a Th1-type T cell. Accordingly, the process may include the steps of: (a) contacting a first sample of mononuclear cells

derived from a patient with OKT3 at or below 37°C to produce a T3CS; (b) removing the T3CS from the first sample; (c) determining the concentration of IFN-γ, IL-4, and IL-10 in the T3CS, and if required, adjusting the concentration of IFN-γ to at least 500 pg/ml, the concentration of IL-4 to less than or equal to 20 pg/ml, and the concentration of IL-10 to less than or equal to 20 pg/ml; (d) providing a second sample of mononuclear cells derived from the same patient; (e) inactivating suppressor cells in the second sample; and, (f) contacting the second sample with T3CS and for a period of time sufficient to activate the second sample in vitro to yield a population of immunoreactive cells.

The process of the invention may also include a 15 step in which the state of activation of the immunoreactive cells is evaluated prior to administration of the cells to the patient. The cells may be evaluated phenotypically or functionally, i.e., by measuring the expression of a cell surface marker indicative of cell 20 activation and/or differentiation or by measuring cell proliferation in response to an additional immune stimulus, e.g, antigen or phorbol myristate acetate (PMA). For example, following cell processing according to the invention, the number of CD25-positive cells in 25 the population of immunoreactive T cells may be measured and the cells discarded if the number is less than 10% of the total number of T cells in the population. Preferably the number of CD25⁺ cells is at least 20% of the total number of T cells in the population. 30 level of activation of the immunoreactive cells may also be evaluated by measuring proliferation in response to further stimulation by immune stimulants. The cells are discarded if the level of proliferation, e.g., the amount of ³H-thymidine incorporated into cellular DNA, is less 35 than twice the level of proliferation of a sample of

unprocessed mononuclear cells. The level of proliferation of the immunoreactive cells is at least twice that of unprocessed cells and preferably is 5-fold greater than that of unprocessed patient-derived cells.

5 If necessary, the level of activation of the immunoreactive cells may be adjusted, e.g., to a higher level of activation, by altering the length of the patient-derive mononuclear cell-T3CS co-culture period to generate the immunoreactive cells, or alternatively, by adding fresh T3CS.

The mixture of nonspecific autologous lymphocyte activators and OKT3, i.e., T3CS, preferably contains: interleukin-1-alpha (IL-1 α), interleukin-1-beta (IL-1 β), interleukin-6 (IL-6), interleukin-8 (IL-8), TNFa, tumor 15 necrosis factor-beta (TNF β), interferon-gamma (IFN γ), granulocyte macrophage-colony stimulating factor (GM-CSF), monocyte/macrophage colony stimulating factor (M-CSF) and OKT3. In preferred embodiments, T3CS contains 12.7 ng/ml (± 10-40%) OKT3 in addition to autologous 20 cytokines in the following amounts (± 10%-40%): IL-1a (105 pg/ml), IL-1 β (1433 pg/ml), IL-6 (808 pg/ml), IL-8 (213 ng/ml), TNF α (570 pg/ml), TNF β (171 pg/ml), IFN γ (14350 pg/ml), M-CSF (1193 pg/ml), and GM-CSF (840pg/ml). IL-2, interleukin-3, IL-4, interleukin-7, IL-10, 25 IL-12, T cell growth factor-beta (TGF β), and granulocytecolony stimulating factor may each be present in T3CS at a concentration of less than 20 pg/ml. Preferably, these cytokines are present at a concentration of less than 5 pg/ml. At least a 20% increase in the number of CD25+ 30 T cells in the first sample of patient-derived mononuclear cells following a T3CS-generation culture compared to a sample of uncultured mononuclear cells is predictive of sufficient production of autologous cytokines.

To avoid differentiation of the patient-derived mononuclear cells into lymphokine-activated killer (LAK) cells, the concentration of IL-2 must be less than 100 units/ml. The concentration of IL-2 in T3CS is preferably less than 50 units/ml, and more preferably in the range of 10-20 units/ml, and most preferably in the range of 1-5 ng/ml (1 unit of IL-2 is approximately equal to 250 pg/ml).

To generate optimal levels of antigen-independent nonspecific lymphocyte activators, production of T3CS is preferably carried out at a temperature greater than 29°C but less than 37°C, e.g., 35°C, for a period of 2 days. Co-culture of patient-derived mononuclear cells with T3CS to generate the immunoreactive cells may also be carried out at sub-physiologic temperature, e.g., a temperature greater than 29°C but less than 37°C.

Following the incubation of cells with T3CS, the cells may be removed from the T3CS and contacted with IL-2, preferably in an amount which is sufficient to bind to 20 at least 25% of the IL-2 receptors on the surface of the immunoreactive cells; more preferably, the amount of IL-2 is sufficient to saturate the IL-2 receptors on the surface of the immunoreactive cells. Contacting the immunoreactive cells with IL-2 is preferably done at 4°C, e.g., during storage or delivery of the cells prior to administration to the patient.

The invention also features a process of producing a population of antigen-specific polyclonal T cells by (a) contacting a first sample of mononuclear cells

30 derived from a patient with OKT3 at or below 37°C to produce a T3CS; (b) removing the T3CS from the first sample; (c) determining the concentration of OKT3 in the T3CS, and if required, supplementing the T3CS with additional OKT3 to achieve a concentration of at least

35 0.1 ng/ml; (d) providing a second sample of mononuclear

cells derived from the same patient; (e) contacting the second sample with T3CS and an antigen for a period of time, e.g., 1-30 days, sufficient to activate the second sample in vitro to yield a population of antigen-specific polyclonal T cells. The concentration of OKT3 in the T3CS is preferably 0.1-1 ng/ml. To achieve the desired concentration, T3CS may be supplemented with OKT3, or OKT3 may be removed from T3CS using methods known in the art, such as chromatography, antibody-mediated depletion, or filtration. The antigen may be in the form of a natural or synthetic peptide, cell extract, a purified antigen, or a recombinantly expressed antigen and may be a tumor antigen, bacterial antigen, viral antigen, or autoantigen.

15 In another aspect, the invention provides an immunoreactive mononuclear cell produced by the inventive process. The cell is preferably a T cell, more preferably a Th1-type T cell. The T cell preferably expresses at least 10%, more preferably at least 75%, and 20 most preferably at least 100% more cell-surface CD25 than an unprocessed mononuclear cell, e.g., a mononuclear cell in a resting state. The cell of the invention is preferably in a primed state, i.e., the cell proliferates at a rate that is at least twice that of an unprocessed 25 T cell when contacted with an immune stimulant. invention also includes a mixture of immunoreactive cells produced by the inventive process at least 75% of which are T lymphocytes, e.g., Thi-type T cells. The cells of the invention can be used to treat any condition 30 characterized by sub-optimal immune responsiveness.

Another aspect of the invention features a process for producing a mixture of autologous nonspecific lymphocyte activators by collecting mononuclear cells from the blood of a patient afflicted with cancer or an infectious disease, inactivating suppressor T cells in

- 12 -

the sample of mononuclear cells, and contacting the mononuclear cells with a compound that binds to the T cell receptor or the T cell receptor-associated CD3 molecule at or below 37°C, e.g., 29-36°C, e.g., 35°C for 5 2 days. The T cell receptor-binding compound may bind to the alpha chain or beta chain of a T cell receptor (or alternatively to the gamma or delta chain). CD3-binding compound is preferably soluble OKT3, but may be any ligand that binds to the CD3 molecule on the 10 surface of the cell. The cells may be contacted with the CD3-binding compound in the absence or in the presence of an antigen; the CD3-binding compound may be removed from the cell culture supernatant following the production of the mixture of autologous cytokines. In preferred 15 embodiments, the mixture contains autologous cytokines in the following amounts (\pm 10%-40%): IL-1 α (105 pg/ml), $IL-1\beta$ (1433 pg/ml), IL-6 (808 pg/ml), IL-8 (213 ng/ml), TNF α (570 pg/ml), TNF β (171 pg/ml), IFN γ (14350 pg/ml), M-CSF (1193 pg/ml), and GM-CSF (840 pg/ml), either in the 20 presence or absence of 12.7 ng/ml (± 10-40%) OKT3.

The invention also includes the immunoreactive cells of the invention together with a pharmaceutically acceptable carrier or diluent for patient administration.

In yet another aspect, the invention features a
25 method of treating a tumor or a viral pathogen in a
patient by administering to the patient the
immunoreactive cells of the invention. A suppressor cell
inhibiting compound, e.g., cimetidine, indomethacin, or
both, may be concurrently administered to the patient.
30 The method may be used to treat any type of cancer
including both solid tumors and hematologic tumors, e.g.,
renal cell carcinoma, breast carcinoma, prostate
carcinoma, colo-rectal carcinoma, pancreatic carcinoma,
ovarian carcinoma, melanoma and non-small cell carcinoma
of the lung as well as leukemias and lymphomas. The

- 13 -

methods and compositions of the invention represent a promising approach to tumors not treatable by conventional forms of therapy such as chemotherapy, radiation therapy, or surgery.

Mononuclear cells taken from a patient afflicted with a complex chronic viral disease may also be processed according to the invention to yield immunoreactive cells which can then be returned to the patient to augment the patient's immune response to the pathogen. Patients infected with pathogenic viruses, e.g., hepatitis B virus, hepatitis C virus, recurrent herpesvirus (herpes simplex virus, varicella zoster virus, cytomegalovirus), papilloma virus, Epstein Barr viurs and HIV (HIV-1 and HIV-2), may be treated in this manner.

Other features and advantages of the invention will be apparent from the following detailed description, and from the claims.

Brief Description

The figures will first be briefly described.

Fig. 1A is a flow cytometric histogram showing the level of cell surface CD25 on T lymphocytes cultured in culture media alone.

Fig. 1B is a flow cytometric histogram showing the 25 level of cell surface CD25 on T lymphocytes activated by OKT3 alone.

Fig. 1C is a flow cytometric histogram showing the level of cell surface CD25 on T lymphocytes activated with 25% (vol/vol) T3CS (patient #1).

Fig. 1D is a flow cytometric histogram showing the level of cell surface CD25 on T lymphocytes activated with 25% (vol/vol) T3CS (patient #2).

- 14 -

Fig. 1E is a flow cytometric histogram showing the level of cell surface CD25 on T lymphocytes activated with 25% T3CS (patient #3).

Fig. 2A is a bar graph showing cell surface CD25

(IL-2 receptor) expression of peripheral blood
mononuclear cells (PBMC) obtained from 7 metastatic renal
cell carcinoma (mRCC) patients following a 5-day culture
with either 25% T3CS, OKT3-depleted T3CS, or OKT3 alone.

Fig. 2B is a bar graph showing PMA-induced 10 proliferation of PBMC following a 5-day culture with 25% T3CS, OKT3-depleted T3CS, or OKT3 alone.

Fig. 3A is a line graph showing the enhanced capacity of chimpanzee EVA cells (EVA #2) to proliferate upon stimulation with PMA compared to unprocessed chimpanzee-derived mononuclear cells (PBMCs).

Fig. 3B is a line graph showing the enhanced capacity of chimpanzee EVA cells (EVA #3) to proliferate upon stimulation with PMA compared to unprocessed chimpanzee-derived mononuclear cells (PBMCs).

Fig. 4 is a line graph showing the enhanced capacity of EVA cells to proliferate and produce cytokines (IFN γ) upon stimulation with PMA compared to unprocessed mononuclear cells (PBMC).

20

Fig. 5 is a line graph showing the enhanced
25 capacity of EVA cells to proliferate upon stimulation
with recombinant IL-2 compared to unprocessed mononuclear
cells (PBMC).

Fig. 6 is a line graph showing the enhancement of PBMC proliferation in response to a recall antigen

30 (influenza) due to the addition of irradiated EVA cells as helper cells.

Fig. 7 is a bar graph showing the enhanced cytolytic function of EVA cells compared to unprocessed mononuclear cells (PBMC).

Processing of peripheral blood mononuclear cells

Mononuclear cells to be processed according to the invention can be obtained from patients, e.g., those afflicted with a malignant tumor or an infectious disease 5 such as hepatitis B. Peripheral blood or a mononuclear cell-enriched population of cells (obtained using known methods, e.g., apheresis) is taken from a patient, and a portion of the sample is mixed with an anticoagulant, e.g., heparin, sodium citrate, ethylenediaminetetraacetic 10 acid, sodium oxalate. The blood-anticoagulant mixture then is diluted in a physiologically acceptable solution such as sodium chloride or phosphate buffered solution. Mononuclear cells are recovered by layering the bloodanticoagulant composition onto a centrifugation 15 separation medium such as Ficoll-Hypaque (Pharmacia Corporation) or Lymphocyte Separation Medium (Litton Bionetics Corporation). The layered mixture then is centrifuged, and the interface containing the mononuclear cells is collected and washed.

The suppressor cells in the mononuclear cell 20 population may be functionally inactivated by contacting the mononuclear cells with an agent that has a specific affinity for or effect upon suppressor cells. particularly suitable composition for inactivating 25 suppressor cells is an H2 receptor antagonist, such as cimetidine; a suitable composition for inactivating the suppressor activity of monocytes is indomethacin. Following suppressor cell inactivation, the mononuclear cells are suspended in a culture medium containing a 30 mitogenic compound which binds to the T cell receptor or the T cell receptor-associated CD3 molecule, e.g., a CD3binding compound, e.g., OKT3, to produce T3CS. Cimetidine may be used in the inventive process to inactivate suppressor cells. The addition of cimetidine 35 to the medium when PBMC are cultured in the T3CS resulted

in increased activation of T cells as measured by enhanced proliferative responses of immunoreactive cells upon further stimulation with PMA (data not shown).

In each sample, the concentration of mononuclear cells can be in the range of about 0.5 - 5.0 x 10⁶ cells/ml, preferably 1-2 x 10⁶ cells/ml. Although any standard tissue culture medium can be utilized in the process of this invention, the cells are preferably cultured under serum-free conditions at 37°C using a standard tissue culture medium, e.g., AIM V medium available from Gibco-BRL, Grand Island, NY.

Mononuclear cells are generally cultured with OKT3 for a period of 3 days to generate T3CS. The T3CS may be used immediately or stored frozen and then thawed for 15 use. The concentration of cytokines, e.g., TNFα, or OKT3 concentration in the T3CS may be measured by any conventional means such as radioimmunoassay or enzymelinked immunosorbent assay (ELISA) using antibodies specific for those components.

20 EXAMPLE 1: Characterization of T3CS Generation of T3CS

Peripheral blood mononuclear cells (PBMC) were obtained from patients with mRCC by leukopheresis and fractionated using ficoll density separation. The cells from the mononuclear fraction were cultured at 1 x 10⁶/ml in AIM V medium (Gibco-BRL, Grand Island, NY) with 25 ng/ml OKT3 (Orthoclone OKT3; Ortho Pharmaceutical Corporation, Raritan, NJ) for 3 days in Lifecell® bags. To inhibit suppressor cell activity, 50 µM cimetidine (Tagamet®; Smith Kline Beecham Pharmaceutical, Cidra, PA) and 10 nM indomethacin (Indocin®; Merck Sharp & Dohme, West Point, PA) were also added to the culture medium. At the end of the culture, the culture bags were centrifuged at 1100 × g for 20 min at room temperature,

and the supernatants were collected, aliquoted, and stored at -70°C.

Composition of T3CS

The composition of the T3CS was determined by 5 ELISA analysis using Quantikine kits from R&D Systems (Minneapolis, MN) for IL-1 α , IL-2, 3, 4, 6, 7, 8, TNF α , TNF β , GM-CSF, TGF β and G-CSF, IFN γ kits from Endogen (Boston, MA) and Gibco (Grand Island, NY), IL-10 kits from Biosource International (Camarillo, CA), and IL-1 β 10 kits from Cistron Biotechnology (Pine Brook, NJ). was determined by a bioassay (phytohemagglutinin (PHA) blast proliferation). The amount of OKT3 present in T3CS was determined by ELISA using the following reagents purchased from Vector Laboratories, Inc. (Burlingame, 15 CA): Horse anti-mouse IgG to capture the OKT3 mouse mAb, biotinylated horse anti-mouse IgG to detect the captured mab and ABC reagent consisting of avidin-conjugated horseradish peroxidase to amplify the signal. Ortho-phenylenediamine dihydrochloride (Sigma Chemical 20 Co., St. Louis, MO) was used as a substrate.

Polyclonal Activation of PBMC

under the same conditions as used in processing patients-derived cells except on a smaller scale (5 ml).

Excess PBMC obtained from mRCC patients were isolated by ficoll density gradient and were cultured at 2 × 10⁶/ml for 5 days in complete AIM V medium (containing 50 M cimetidine and 10 nM indomethacin) with either 25% (vol/vol) T3CS, or various concentrations of OKT3. Cells cultured with medium alone served as control. After incubation, the cells were washed and resuspended at a concentration of 10 × 10⁶/ml in "infusion medium" (1% human serum albumin and 0.5% dextrose in Lactated Ringers

solution). The cells were then stored overnight at 4°C prior to analysis (to simulate overnight storage for final QA/QC and shipping to clinical sites).

Phenotypic Analysis of PBMC and EVA Cells

Cell phenotypes were determined by flow cytometry. 5 Immunofluorescent staining was carried out using mAb specific for cell surface antigens, e.g., Coulterclones T3-RD1/CD3-RD1, T4-RD1, T8-RD1 and -FITC, IL-2R-FITC, I3-FITC and 2H4-RD1 (Coulter Corporation, Hialeah, FL), 10 and Dako UCHL1-FITC (Dako corporation, Carpintera, CA). Appropriate isotype-matched labels were used as negative Cells were resuspended at 2 x 10⁶/ml in AIM V controls. medium and aliquoted to 100 μ l per tube. Two to four μ ls of labeled mAbs were added to the cells according to the 15 manufacturers' recommendation. The cells were then incubated with the Abs at 4°C for 30 min, washed once with 500 μ l of cold PBS, and resuspended at a concentration of 4 x 10⁵ cells/ml in PBS for immediate analysis using a Coulter Epics Profile II flow cytometer.

20 PMA Assay

Proliferative responses to PMA were used to measure the degree of cell activation. The PBMC or EVA cells were resuspended to 1 x 10⁶ cells/ml in AIM V medium with or without 1 ng/ml PMA, and cultured in 25 triplicate in 96-well flat-bottom plates (Costar, Cambridge, MA) at 37°C with 5% CO₂ for 48 hours. The cells were pulsed with ³H-thymidine (1 µCi /well) for the last 6 hours of culture, and harvested onto filtermats. The amount of radioactivity incorporated into the cells was determined by liquid scintillation counting.

Depletion of OKT3 or Cytokines From T3CS

T3CS was divided into 6 ml aliquots, and incubated with the appropriate neutralizing antibody (Goat anti-human, from R&D Systems) at 37°C for one hour. amount of antibody used for depletion was determined by 5 the level of the particular cytokine known to be present and the antibody activity required for neutralization as specified by the manufacturer. Specifically, 125 μ l of antibody was added for the depletion of IL-1 β or TNF α , 20 μ l for IL-6, 60 μ l for IFN, 100 μ l for GM-CSF, and 250 μ l 10 for IL-8. After incubation, 1 ml of pre-washed magnetic beads (Advanced Magnetics, Cambridge, MA) conjugated with either rabbit-anti-goat IgG (for cytokine depletion) or goat-anti-Mouse IgG (for OKT3 depletion) was added to the samples (2 ml of beads were used for multiple-cytokine 15 depletion) and incubated with the cells at room temperature for 20 min. The cytokine-bound beads were then removed by placing the tube in a magnetic holder for 5 min. The cytokine-depleted T3CS was carefully transferred to a fresh labeled tube. The depletions were 20 monitored by ELISA analysis.

Cytokine Composition of T3CS

As described above, T3CS was generated by culturing PBMC isolated from mRCC patients with 25 ng/ml OKT3 for 3 days. Aliquots of randomly selected mRCC patient T3CS (N=33-42) were analyzed for the presence of 17 different cytokines and OKT3 by ELISA. As shown in Table 1, T3CS contained a mixture of monokines and lymphokines including IL-1α, IL-1β, IL-6, IL-8, TNFα and β, IFNγ, M-CSF, and GM-CSF. Several other cytokines, including IL-2, 3, 4, 7, 10, 12, TGFβ, and G-CSF, were either undetectable or detected at very low levels in all samples tested. The cytokine profile indicated that the cells involved in the PBMC activation process were predominantly monocytes and Th1-type T cells (Seder et

al., 1994, Ann. Rev. Immunol. 12:635-673) or T cells that had differentiated into Th1-type cytokine producers during the cultures. Although kinetic studies indicated that significant amounts of IL-2 (15-433 pg/ml) were 5 present in all samples of T3CS analyzed during the first 24 hours of the culture, the level of IL-2 dropped sharply by day 2 (0-132 pg/ml). IL-2 levels became undetectable by day 3 when the culture supernatant was harvested. The drop in IL-2 is likely to be due to 10 active consumption of the cytokine during cell culture. The mean OKT3 concentration in T3CS was 12.7 ng/ml. contains significantly lower levels of IL-2 than conventional conditioned media, e.g., PHA-generated cell supernatants, Concanavalin-A-generated cell supernatants, 15 or mixed lymphocyte culture supernatants. Cytokine levels in T3CS from the majority patients fell into acceptable ranges, i.e., sufficient amounts of the key cytokines were present to produce immunoreactive cells in subsequent PBMC cultures. T3CS from a small number of 20 patients contained low levels of certain cytokines. As discussed above, the concentration of cytokines may be augmented, if necessary, by carrying out the T3CS generating step for a longer period of time or alternatively, by adding recombinant or nonrecombinant 25 cytokines.

- 21 -

TABLE 1

COMPOSI	TION OF CONDITION	ONED MEDIUM (1	r3C8)
		Mean	N
Cytokines:			
IL-1α	(pg/ml)	105	33
IL-1β	(pg/ml)	1433	36
IL-6	(pg/ml)	808	39
IL-8	(ng/ml)	213	39
TNFα	(pg/ml)	570	39
TNFβ	(pg/ml)	171	39
INFγ	(pg/ml)	14350	40
GM-CSF	(pg/ml)	840	39
M-CSF	(pg/ml)	1193	39
октз	(ng/ml)	12.7	42

IL-2, 3, 4, 7, 10, 12, $TGF\beta$ and G-CSF are below detectable levels (<3-8 pg/ml) in >90% samples tested.

Accessory Cell Requirement for Production of Immunoreactive Cells

10

According to the invention, the concentration of patient-derived monocytes in a sample of patient-derived mononuclear cells, is preferably 0.1-50%, more preferably 1-30%, more preferably 5-15%, and most preferably 10%.

To evaluate the role of monocytes in cultures of patient-derived mononuclear cells, monocytes were depleted from patient-derived mononuclear cells using the well-known methods of adherence, e.g., to plastic plates, or incubation with L-leucine methyl ester. Incubation of a monocyte-depleted sample of patient-derived mononuclear cells with T3CS failed to yield the immunoreactive cells of the invention. In addition, blocking of the Fc receptor on the surface of Fc-receptor positive accessory

cells, e.g., by adding an excess of soluble human polyclonal IgG to the patient-derived mononuclear cell-T3CS co-culture, inhibited the T3CS-catalyzed generation of immunoreactive cells by 93%. These data indicate that Fc-receptor positive accessory cells, i.e., monocytes, granulocytes or dendritic cells, are required to generate immunoreactive cells using the methods of the invention when OKT3 is used in solution phase.

If a sample of patient-derived mononuclear cells has a sub-optimal concentration of accessory cells, the sample of patient-derived mononuclear cells may be enriched for monocytes, granulocytes, or dendritic cells.

Role of Soluble OKT3

10

15

20

30

The OKT3 used in the process of the invention is preferably in solution phase rather than solid phase. One advantage of using soluble OKT3 in the inventive process is that soluble OKT3 mediates a more physiological interaction between Fc-receptor-bearing accessory cells, e.g., monocytes, and T cells. By forming a bridge between these cells, a full complement of costimulatory signals is initiated, thus minimizing the potential for the generation of incompletely (anergic) or aberrantly (apoptotic) activated T cells.

Enhancement of OKT3-Induced T Cell Activation by Autologous Cytokines

The role of cytokines in the induction of T cell activation was compared to that of OKT3 alone. PBMC from seven mRCC patients were cultured for 5 days in AIM-V medium containing either 2.5 ng/ml OKT3 alone or 25% (vol/vol) T3CS that was first depleted of OKT3 and then reconstituted with 2.5 ng/ml of fresh antibody. This depletion/reconstitution approach was undertaken to assure that the amount of biologically active OKT3 in the

5

10

15

20

25

30

T3CS would be identical to the OKT3 control. To avoid the bias from a single aliquot of T3CS, a pool made from three different mRCC patients was used to stimulate the PBMC from all seven patients.

The level of T cell activation was determined by the expression of cell surface IL-2 receptor (CD25) as measured by flow cytometry. Figs. 1A-1E are representative immunofluorescence histograms that demonstrate that CD25 expression on T cells stimulated with the autologous cytokine-OKT3 mixture was significantly higher than that on T cells stimulated with OKT3 alone.

Primed Activation State of Immunoreactive Cells

The degree of T cell activation was also assessed functionally by measurement of the proliferative response of EVA cells upon further stimulation by PMA, a protein kinase C activator.

Fig. 4 shows the results of an experiment in which immunoreactive cells were contacted with an immune stimulant, e.g., various concentrations of PMA. absence of further stimulation by PMA, immunoreactive cells displayed very little spontaneous proliferation or cytokine secretion when cultured at 37°C in AIM V medium alone. However, when contacted with 1-10 ng/ml of PMA, immunoreactive cells (but not unprocessed PBMC) displayed an enhanced capacity to proliferate and produce IFNy (as well as other Th1-type cytokines such as $TNF\alpha$, $TNF\beta$, and GM-CSF (data not shown)). These results are in contrast to those obtained with mononuclear cells activated by a conventional stimulant, such as PHA, which possess high levels of spontaneous proliferation and cytokine production following an initial 5-day activation culture (data not shown).

- 24 -

This enhanced proliferation and cytokine secretion upon contact with an immune stimulant and concomitant absence of or low levels of spontaneous proliferation and cytokine secretion indicates that the immunoreactive cells of the invention are in a primed state of activation.

Similarly, an increase in proliferation was observed when EVA cells were contacted with another immune stimulant, i.e., IL-2 (Fig. 5). These data confirm the primed nature of the immunoreactive cells and suggest that the immunoreactive cells of the invention could be effectively co-administered with a low dose of IL-2.

Enhanced Effector Function of EVA Cells: Cytolytic T Cells

10

30

The ability of the T3CS to support the generation of cytótoxic EVA cells was evaluated. As shown in Fig. 7, cytotoxicity was determined using a conventional chromium-51 release assay with the K562 leukemia cell line and allogeneic human renal carcinoma cell line 769P as targets. As shown in Fig. 7, EVA cells, i.e., T cells which have been polyclonally-activated independent of tumor-associated antigens according to the invention, possessed a greatly enhanced cytotoxicity toward both tumor lines compared to unprocessed PBMC. These in vitro results suggest that EVA cells are capable of directly killing tumor cells in vivo.

Enhanced Effector Function of EVA Cells: Helper T Cells

In addition to enhanced cytolytic function, the immunoreactive cells of the invention were found to have enhanced helper function. The ability of T3CS to support the generation of EVA cells with helper cell function was determined by measurement of the ability of irradiated

PCT/US95/00193 WO 95/20649

EVA cells to enhance the proliferative response of unprocessed mononuclear cells upon stimulation by a recall antigen, e.g., an influenza antigen. As shown in Fig. 6, EVA cells added to cultures at low levels (5-10%) provided helper signals to unprocessed patient-derived mononuclear cells (presumably through the secretion of Th1-type cytokines) resulting in a significant increase in proliferation. These in vitro results suggest that EVA cells are capable of amplifying and broadening their effects in vivo through the production of cytokines.

10

15

20

25

The ability of EVA cells to proliferate and to produce a variety of cytokines (IL-2, GM-CSF, IFN γ , TNF α) in vitro in response to further stimulation by such agents as PMA and IL-2, as well as to lyse tumor cell targets, is greatly enhanced compared to the PBMC from which they were derived. The lowered activation threshold of the EVA cells exhibited in vitro suggests that once they are reinfused into patients, they are likely to demonstrate enhanced responsiveness to immunological signals, such as weakly immunogenic tumor antigens which normally are non-stimulatory to unprocessed cells.

Phenotypic Characterization of Immunoreactive Cells

Following short term (5-day) culture of patientderived mononuclear cells in the autologous cytokine mixture/OKT3-containing conditioned medium, the resulting EVA cells were analyzed for the expression of cell surface antigens. EVA cells expressed enhanced levels of a variety of activation and/or differentiation markers on their cell surface including cytokine receptors, e.g., 30 CD25, major histocompatibility complex (MHC) antigens, e.g., MHC class II, adhesion molecules/homing receptors, e.g., CD44/Leu8, costimulatory molecules, e.g., CD28, and markers of primed or memory T cells, e.g., CD45RO, as shown in Table 5.

TABLE 5

EVA CELL PRODUCT IDENTITY

	%(+) LYM	PHOCYTES
	PBMC	EVA CELLS
CD3(+) - T CELLS	75%	84%
CD4(+) - HELPER T CELLS	48%	60%
CD8(+) - CYTOTOXIC T CELLS	20%	25%
CD3(+)/CD25 ⁽⁺⁾ - ACTIVATED T CELLS (IL2 RECEPTOR-EARLY STAGE MARKER)	3%	44%
CD3(+)/CLASS II(+) - ACTIVATED T CELLS (MHC-LATE STAGE MARKER)	13%	42% ·.
CD3(+)/CD45RO(+) - T CELLS (PRIMED "MEMORY" CELLS)	49%	69%
CD45RO(+)/CD25 ⁽⁺⁾ - LYMPHOCYTES	3%	36%
CD3(+)/LEU 8(+) - T CELLS (LYMPH NODE HOMING RECEPTOR)	31%	59%
CD3(+)/CD28(+) - T CELLS (B7/BB1 RECEPTOR - COSTIMULATORY SIGNAL TRANSDUCTION)	0.9*	2.3*
CD3(+)/CD44(+) - T CELLS (EARLY ACTIVATION MARKER - CELL ADHESION MOLECULE)	6.0*	8.3*

*: Mean Channel Fluorescence Intensity N = 19 mRCC Patient PBMC & EVA Cell Products

Generation of Immunoreactive Cells: Synergistic Effects of Autolog Cytokines Plus OKT3

OKT3 was depleted from the T3CS in order to determine the relative effects of the anti-CD3 monoclonal

antibody and the autologous cytokines on the induction of CD25 expression on the surface of T lymphocytes (Fig. 2A) and the proliferation of EVA cells upon further stimulation by PMA (Fig. 2B). Following depletion of OKT3 from the T3CS, there was little or no measurable activation of T cells, i.e., the autologous cytokines were not capable of stimulating resting PBMC in the absence of OKT3. However, when mRCC patient PBMC were cultured with complete T3CS, the cytokines functioned in synergy with OKT3 resulting in large increases in T cell activation relative to the levels achieved with OKT3 alone. Taken together, these results indicate that OKT3 catalyzes the generation of immunoreactive cells, while the autologous cytokines function as costimulants to optimize the overall T cell activation process.

Immunoreactive Cell Generation Utilizing Various Recombinant Cytokine-OKT3 Mixtures

10

15

20

30

Reconstitution experiments were carried out to identify the cytokines in the T3CS have the greatest effects on T cell activation. The roles of 5 major cytokines (IL-1 β , IL-6, TNF α , IFN- γ , GM-CSF) were analyzed utilizing various combinations of recombinant cytokines together with OKT3. Table 2 shows the results of a series of experiments in which PBMC from 7 mRCC cancer patients were cultured in an allogeneic T3CS, OKT3 alone, or various OKT3-multiple recombinant cytokine mixtures. Five out of the six cytokines analyzed were found to enhance the level of stimulation of T cells compared to the level of stimulation seen with OKT3 alone, as measured by cell surface expression of CD25 on T cells. This result is consistent with the observation that the cytokines in T3CS act synergistically with the OKT3 catalyst in the activation process. Upon analyzing the data from individual mRCC patients, however, no

- 28 -

single cytokine could be identified as the dominant costimulant across all cultures. In addition, for the majority of patient PBMC cultures, optimal T cell activation was achieved only with complete conditioned medium (T3CS). Taken together, these results indicate that a single recombinant cytokine, or even a mixture of several recombinant cytokines, is not able to provide the full complement of costimulatory signals required for optimal polyclonal activation of T cells. In contrast, use of a multiple autologous cytokine-enriched conditioned medium, i.e., T3CS, consistently resulted in the generation of optimally immunoreactive cells.

10

TABLE 2

Comparison of Conditioned Medium with OKT3 Plus Recombinant Cytokines on T Cell Activation in mRCC patients

		% CD25	% CD25 Positive Cells	Cells				
Stimulants	,	-		πВ	mRCC Patient			
	· Average	# 1	#2	#3	# 4	#2	9#	L#
Medium alone	7	6.4	5.4	5.3	5.4	4.9	5.5	14.0
OKT3 2.4 ng/ml	38	18.8	24.7	54.7	28.6	31.7	64.1	46.3
OKT3 + IL1* + TNFa*	4 1	40.4	37.3	58.1	31.3	36.3	71.8	48.2
OKT3 + IL1 + IFNg*	44	32.3	37.1	59.7	16.4	39.1	71.9	51.9
OKT3 + IL1 + GM-CSF*	42	26.9	38.0	55.9	22.8	38.1	6.69	43.2
OKT3 + IL1 + IL6*	48	43.5	43.2	55.3	26.5	37.5	70.8	54.9
OKT3 + IL1 + IFNg + GM-CSF	45	33.5	44.4	58.4	21.6	34.8	68.7	41.9
OKT3 + IL1 + TNFa + IL6	43	20.4	35.0	69.3	24.1	33.7	70.1	44.2
OKT3 + IL1 + TNFa + GM-CSF	43	26.1	40.7	66.3	21.8	33.3	67.5	42.2
OKT3 + IL1 + IFNg + IL6	43	26.5	38.1	65.0	21.8	38.2	68.5	47.3
OKT3 + IL1 + TNFa + GM-CSF + IL-6	41	27.6	39.7	61.2	18.8	39.7	73.1	40.8
OKT3 + IL1 + TNFa + IFNg + IL-6	41	40.1	34.6	55.8	16.3	37.5	8.69	42.1
OKT3 + IL1 + TNFa + IFNg + IL-6 + GM-CSF	39	26.1	35.2	39.6	16.0	37.4	72.8	43.9
25% Çöndiliöned Medium**	5.4	53.3	56.4	69.4	57.3	33.2	74.8	33.8

PBMC from seven individual mRCC patients were cultured for 5 days with OKT3 alone, or OKT3 plus various combinations of cytokines, or with 25% of an allogeneic conditioned medium. The cells were stored at 4°C overnight and then analyzed for CD25 expression by flow cytometry.

**25% condilioned medium (WP35) contains: OKT3 2.4 ng/ml, IL1b 650pg/ml, IL6 400pg/ml, TNFa 131 pg/ml, IFNg 1341pg/ml, GMCSF 57pg/ml. *Concentration of recombinant cytokines: IL-1b=250 pg/ml, IL-6=100 pg/ml, TNFa=100pg/ml, IFNg=250pg/ml, GM-CSF=250pg/ml.

TNFa Is an Important Co-stimulant for T Cell Activation

Table 3 shows the results from a series of experiments in which the contributions of various cytokines present in the T3CS to the overall T cell activation process was analyzed. The level of induction of CD25 on PBMC cultured in T3CS was compared to that achieved when cytokines were selectively depleted from aliquots of the same T3CS preparation.

Despite quantitative differences between the various RCC patient PBMC tested, depletion of TNF α from the T3CS consistently resulted in a reduction of CD25 expression on the EVA cells generated from all seven cancer patients, indicating that TNF α serves a key costimulatory function in enhancing OKT3-catalyzed polyclonal T cell activation.

10

15

20

25

30

The effects of the depletion of other cytokines from the T3CS were highly variable suggesting that the requirement for other costimulatory signals for optimal T cell activation varies from patient to patient, and perhaps even from culture to culture dependent upon the quality of the PBMC isolated at various points in time from different patients at different times along their disease course and upon the expression of cytokine receptors on freshly isolated mononuclear cells.

These results confirm that a single recombinant cytokine, or even a mixture of several recombinant cytokines, is not able to provide the necessary costimulatory signals for optimal polyclonal activation of T cells, and the complete spectrum of autologous cytokines in T3CS is necessary to generate the immunoreactive cells of the invention.

TABLE 3

Effect of Depletion of Various Cytokines from Conditioned Medium on T Cell Activation

,		IL-2	Receptor Ex	IL-2 Receptor Expression (% of Control)	s of Control).	
Patient #	Complete	***************************************		Cytokine Depleted	Depleted	eted	***************************************
	Conditioned Media	TNFa	11-1	IL-6	IL-8	GMCSF	FINA
-	100	84	88	110	7.8	96	66
2	100	88	92	115	77	၉	98
က	100	6.2	91.	86	68	8.7	103
4	100	85	95	109	85	122	88
ľ	100	7.1	92	86	66	66	185
.	100	65	105	. 92	80	85	140
7	100	68	83	85	117	117	115
Average	100	7.7	0 6	9.8	8 6	95	118

PBMC from seven different mRCC patients were stimulated with 25% allogeneic conditioned media either complete or depleted with (% CD25 positive cells x mean channel fluorescence intensity) on cells activated by the complete, undepleted conditioned media. one of the cytokines listed. The cells were harvested 5 days later and examined for surface IL-2 receptor (CD25) expression by flow cytometry the next day. The effect of cytokine depletion is expressed as % of control which is the total IL-2R expression

- 32 -

T3CS was utilized as a nonspecific stimulant in the ex vivo generation of polyclonally activated T cells for the treatment of metastatic renal cell carcinoma.

T3CS produced according to the invention contains a wide variety of monokines and lymphokines together with low levels of OKT3. OKT3 functions synergistically with the cytokines in the T3CS. Removal of OKT3 from the T3CS resulted in a substantial decrease of T cell activation demonstrating the role of anti-CD3 monoclonal antibody as a catalyst while the cytokines provide the costimulatory T cell activation signals.

10

15

20

25

30

35

Conditions in which OKT3 and autologous cytokines function synergistically in T cell activation have been achieved in serum-free medium. Serum-free culture is particularly desirable for the generation of EVA cells used in adoptive immunotherapy.

The claimed process employs aliquots of the T3CS . generated in advance to be used as a stimulant in the secondary cultures. This approach ensures that a full complement of costimulatory signals are available at the initiation of each activation culture and therefore minimizes the probability of generating anergic or apoptotic T cells. Furthermore, use of the pre-manufactured and quality-assured autologous cytokine-OKT3 mixture as a stimulant decreases the dependence upon de novo synthesis of cytokines during the early stage of the subsequent T cell activation cultures. Such decreased dependency is especially important under the conditions such as off-site cell processing that require shipping and storage of the cells, and when dealing the with cells from patients of various clinical stages.

PBMC obtained from individual cancer patients, and even PBMC taken from the same patient at different times, respond differently to various cytokines. OKT3-induced

- 33 -

T cell activation for one patient's cells may strongly be affected by a given cytokine while cells from another patient are not affected at all by the same cytokine. This variation is due both to the inherent diversity of human lymphocytes and accessory cells, and to the fluctuation caused by the process prior to the activation such as collection of the cells by apheresis. Overnight storage and shipping may affect accessory cell activity which in turn affects the dependence of the T cells on exogenous cytokine present in the culture. Hence, the methods of the invention require evaluation of various key components of the T3CS, e.g., OKT3 and TNF α concentration, and the concentration of accessory cells, e.g., monocytes. The methods may also require the evaluation of the activation state, i.e., PMA responsiveness, of the processed cells to predict their clinical effectiveness.

10

15

20

25

30

IL-2-Independent Generation of Immunoreactive Cells

Utilizing T3CS as a stimulant for mononuclear cells, it is possible to generate T cells expressing high levels of CD25/IL-2 receptors, independent of the inclusion of high levels of exogenous IL-2 in the culture medium. Although low levels of IL-2 are present in both the T3CS and EVA cultures at early timepoints, there is no detectable (< 6 pg/ml) autologous IL-2 present in the T3CS when the EVA culture is initiated. Consequently, in contrast to LAK and TIL cells which are both cultured in high levels of IL-2, the immunoreactive cells generated by the inventive process are presumably less dependent upon systemic administration of high dose IL-2 for therapeutic efficacy. In addition, the lack of IL-2 in the T3CS allows the production of immunoreactive cells, i.e., multifunctional, polyclonal T cells containing both CD8⁺/cytotoxic T cells and a high percentage of

CD4⁺/helper cells, in contrast to PBMC grown in high dose IL-2 which are highly enriched in CD8⁺/cytotoxic T cells. Thus, the immunoreactive cells of the invention have broader functional capacities than PBMC cultured in IL-2.

If production of T3CS containing low levels, e.g., up to 10 units/ml, of autologous IL-2 is desired to further enhance the generation of immunoreactive cells, IL-2 production may be increased and IL-2 consumption simultaneously decreased by modification of the culture conditions, such as the kinetics and temperature of the T3CS generation process.

Following processing, the immunoreactive cells may be exposed to IL-2 at 4°C such that IL-2 binds to at least 25% of the cell surface IL-2 receptors.

Preferably, the cell surface IL-2 receptors are saturated, i.e., 100% of the IL-2 receptors on the surface of a cell are bound to IL-2, with IL-2 prior to infusion into a patient.

Immunoreactive cells with cell surface bound IL-2 are likely to have an enhanced ability to expand in vivo and a decreased dependence upon helper cell-mediated IL-2 production, an activity which may be lacking or depressed in immunosuppressed cancer patients. In addition, administration of IL-2 to a patient in a cell-bound form avoids the toxicity and other clinical complications often associated with intravenous or subcutaneous co-administration of high dose IL-2 to support cell therapy.

OKT3

5

10

20

25

Intravenous OKT3 administration, utilized for the suppression of transplant graft rejection in humans, induces rapid (within 24 hours) removal of a high percentage of CD3⁺ T cells from the circulation due to the presence of the foreign antibody molecules on the

lymphocyte cell surface (Vigeral et al., 1986, Transplantation 41:730; Chatenoud et al., 1986, J. Immunol. 137:830-838). Excessive amounts of OKT3 bound to an EVA cell surface may accelerate the clearance of the T cells, resulting in decreased therapeutic efficacy. The use of T3CS as a stimulant minimizes this clinical problem because T cell activation can be achieved with a minimal concentration of OKT3 in the presence of autologous cytokines compared to the concentration required in the absence of autologous cytokines.

10

15

20

25

30

The use of the T3CS to stimulate PBMC allows a high level of T cell activation to be catalyzed by a minimal amount of OKT3 (1-4 ng/ml). The inventive methods leave very little detectable surface-bound antibody on the activated T cell product, and therefore decrease the probability of the EVA cells being rapidly removed from the circulation. In addition, the low level of surface-bound OKT3 allows more CD3-TCR complexes on the surface of the activated T cells to remain unoccupied so that key molecules such as MHC Class I and II molecules, tumor antigens, and immunogenic peptides may bind to and stimulate the T-cell receptor following re-infusion of the cells into patients. Stimulation of PBMC with higher concentrations of OKT3 tends to drive the cell culture towards CD8 cell dominance, whereas culture in the autologous cytokine-OKT3 mixture maintains a relatively high CD4/CD8 ratio which is desirable for the therapeutic function of EVA cells. Finally, the low level of cell-bound OKT3 reduces the danger of a patient developing HAMA responses after multiple infusions of EVA cells.

EXAMPLE 2: Effect of Temperature on Cytokine Production Conditions under which patient-derived mononuclear cells are cultured with OKT3 to generate T3CS were

- 36 -

analyzed. For the generation of T3CS, cells were cultured for a period of 1 to 5 days. As shown in table 4, it is possible to manipulate the time course and/or temperature at which T3CS is generated in order to selectively increase or decrease the levels of certain autologous cytokines.

The effect of temperatures ranging from 29°C to 39°C on the production of cytokines was also tested. As shown in Table 4, higher levels of cytokines were detected in cultures grown at temperatures of 37°C or below.

10

17.2 16.2 15.5

711.2 542.2

270.5

365.1

601.7

071.0

272.6 350.2

166.3

1170.0

1258.6

421.5

621.4

709.1

122.2

319.5

707.9

204.3

74.8

381.1

11-2

TNFå

11.6

434.9

TABLE 4

1006.9 525.5 3043.8 633.7 270.8 0.069 1393.7 1578.1 3192.4 563.0 1165.3 210.1 43.4 946.1 944.0 363.1 466.7 108.1 19.7 17.9 9999 AVERAGE OF ALL FIVE PATIENTS 593.6 1634.6 2460.3 403.8 1146.4 1579.6 1179.5 758.8 53.0 1378.8 530.0 409.5 1371.5 311.2 289.5 975.6 0.0 17.0 15.9 12.8 10.0 2222 1435.3 1540.3 1430.0 568.3 813.1 672.5 751.3 1463.9 822.0 533.5 78.9 231.5 3101.5 406.1 432.4 3086.2 633.3 2.0 497.2 21.3 20.9 17.4 12.1 2222 1036.1 478.3 3039.3 5232.0 1658.0 1498.6 1195.2 37.0 31.7 804.5 741.4 4009.2 795.0 304.7 522.2 715.8 397.1 322.2 993.1 420.7 2.5 205.1 455.1 0.00019.2 15.9 11.9 7.6 183.9 469.6 39°C 569.9 276.9 1012.0 1680.3 17.4 149.7 422.8 653,9 626.2 542.0 416.6 401.9 345.6 348.2 949.2 603.6 370.3 535.2 3.5 13.6 19.5 14.2 Temperature/ Timecourse Study Days 2 6 2 2 6 5 5 3 200 205 2 6 8 1861.5 1076.0 1440.0 29°C 831.6 689.3 228.0 2706.6 198.4 597.6 592.2 660.1 233.6 258.5 407.4 17.6 991. 839 453 13.5 10.0 6.01 297 2 2 2 2 2 3 3 3 9999 9999 9999 9999 9999 **용 및 및 및** 문문문문 1581.6. 2134.8 1735.2 1100.0 946.3 441.2 702.9 237.0 351.0 506.5 605.5 341.6 1,532 994 1,527 13.9 11.8 8.8 6.0 1.8 34°C 문문문 9999 9999 물 물 물 물 9 9 9 9 5 5 5 5 물물물물 1819.5 1318.2 1537.6 557.8 1774.2 1486.4 1438.4 18.6 830.4 740.0 497.6 590.6 1.3 3.2 782 1,862 1,747 1,653 13.0 2222 6.3 1031.2 1427.1 1237.2 1064.4 1740.8 521.4 629.4 1879.2 68.9 292.8 416.8 549.3 925 1,916 1,675 1,780 667.2 6.0 9.0 9.4 2222 Days 2 6 2 2 6 2 2 6 6 Cytokine 1,

GMCSF

IL-113

giFN

OKT3

EXAMPLE 3: Production and therapeutic administration of EVA cells

Human peripheral blood cells are collected by apheresis, and the PBMC are isolated using a Ficoll gradient. The PBMC are incubated at 37°C at a concentration of 10⁶ cells/ml in AIM V medium containing cimetidine (5 x 10⁻⁵ mol/l), indomethacin (10⁻⁸ mol/l), and 25 ng/ml OKT3 (Ortho Biotech, Raritan, NJ). The activity of suppressor T cells and of suppressor monocytes is inactivated by culturing the PBMC with cimetidine and indomethacin. After three days, the culture supernatant, i.e., T3CS, is harvested, divided into portions, and frozen at -80°C until required for use in the cell cultures.

5

10

patient, and suspended at approximately 2 x 10⁶ cells/ml in AIM V medium containing 25% v/v T3CS, cimetidine (5 x 10⁻⁵ mol/l) and indomethacin (10⁻⁸ mol/l). The cells are harvested following a five day culture at 37°C in a 20 moist-air incubator containing 5% CO₂, then resuspended at 10⁷ cells/ml. Following overnight storage of the activated cells at 4°C, approximately 10⁹ cells are reinfused intravenously at room temperature to the patient over 30 min.

25 If the concentration of OKT3 or TNFα or other cytokines in the T3CS are below the preferable ranges, the levels of one or more of the cytokines may be increased as follows. One way to increase the amount of a cytokine, e.g., TNFα, is to adjust, e.g., increase, the 30 time of culture before harvesting the T3CS. Alternatively, the concentration of monocytes (which produce TNFα) may be increased in the T3CS-generating culture, or the temperature at which the T3CS-generating step is carried out may be adjusted to optimize TNFα production. Another way to increase the amount of a

10

15

20

25

30

cytokine is to increase the percentage of T3CS added to the second sample of patient-derived mononuclear cells, i.e., by adding greater than 25% T3CS.

The therapeutic methods of the invention which utilize EVA cells may be used to treat mRCC and other types of metastatic cancers, as well as infectious diseases, autoimmune diseases, and immunodeficiency diseases. Clinical outcome may be assessed by such measures as length of patient survival, quality of life measurements, changes in any indicators of medical function such as clinical chemistries, size of tumors, changes in load of virus, bacteria, fungus, or parasite, toxicity of the therapy, or delay in time of recurrence of the disease, or other assessments. Clinical outcome may also be evaluated by monitoring changes in immune status or time to recurrence of a tumor in a cancer patient.

EXAMPLE 4: Treatment of Hepatitis B Infection

Polyclonal Activation of PBMC from Chimpanzees with Chronic Hepatitis B Infection

Viral hepatitis refers to an infection of the liver caused by a small group of hepatotropic viruses that may be differentiated serologically. Of the three serotypes, A, B and C; B and C cause far more serious infections than hepatitis A.

The methods of the invention were used to produce immunoreactive cells from mononuclear cells derived from chimpanzees infected with hepatitis B, a primate model of chronic viral infection. The chimpanzee model was chosen over the woodchuck or mouse model as the most relevant model of the human disease because chimpanzees can be experimentally infected with the human hepatitis B virus.

Mononuclear cells from chimpanzees with chronic hepatitis B were polyclonally activated with OKT3 to

- 40 -

produce a T3CS. A second sample of mononuclear cells was cultured with T3CS to obtain EVA cells. Chimpanzeederived mononuclear cells processed according to the invention yielded immunoreactive cells. As shown in Figs. 3A and 3B, the immunoreactive cells were in a primed state of activation, i.e., they proliferated upon stimulation with PMA whereas unprocessed chimpanzeederived mononuclear cells did not. In addition, the spectrum of phenotypic markers on the chimpanzee EVA cells was found to be similar to that on EVA cells 10 produced from patient-derived mononuclear cells (data not shown). These data indicate that the methods of the invention can be used to generate immunoreactive cells from animals, e.g., patients, with chronic viral 15 infections. These immunoreactive cells can then be reinfused into the patient to augment the patient's immune response to the pathogenic virus.

EXAMPLE 5: In vivo Use of T3CS

30

vivo delivery of cytokines. T3CS which has been depleted of OKT3 may be encapsulated in an appropriate slow release matrix, e.g., a liposome or biocompatible polymer, in order to maintain a critical level of cytokine over time and to minimize the toxicity

25 associated with bolus-delivery of cytokines. In vivo delivery of T3CS may also be targeted to a particular site in the body for optimal effectiveness.

In vivo administration of T3CS offers several advantages over conventional cytokine therapy, i.e., administration of one cytokine. T3CS is a unique mixture of both monokines and lymphokines which initiate and promote activation and/or recruitment of several different immune cell types. Thus, T3CS can induce a full spectrum of immune responses in a patient.

- 41 -

Other embodiments are within the following claims.

- A process of producing a population of immunoreactive cells comprising the steps of:
- (a) contacting a first sample of mononuclear cells derived from a patient with OKT3 at or below 37°C to produce an OKT3-derived culture supernatant (T3CS);
- (b) removing said T3CS from said first
 sample;
- (c) determining the concentration of OKT3 in said T3CS, and if required, supplementing said T3CS with 10 said OKT3 to achieve a concentration of at least 0.1 ng/ml;
 - (d) providing a second sample of mononuclear cells derived from said patient; and
- (e) contacting said second sample with said 15 T3CS and said OKT3 for a period of time sufficient to activate said second sample in vitro to yield a population of immunoreactive cells.
 - 2. The process of claim 1, wherein said OKT3 concentration is in the range of 1-25 ng/ml.
- 20 3. The process of claim 2, wherein said OKT3 concentration is in the range of 10-15 ng/ml.
 - 4. The process of claim 1, wherein step (a) is carried out for 6 hours-7 days.
- 5. The process of claim 4, wherein step (a) is 25 carried out for 3 days.
 - 6. The process of claim 1, wherein said period of time in step (e) is 1-30 days.
 - 7. The process of claim 6, wherein said period of time in step (e) is 5 days.

- 8. The process of claim 1, further comprising after step (d), contacting said second sample with a suppressor cell inhibitory compound.
- 9. The process of claim 8, wherein said suppressor cell inhibitory compound is cimetidine, indomethacin, or both cimetidine and indomethacin.
 - 10. A process of producing a population of immunoreactive cells comprising the steps of:
- (a) providing a first sample of mononuclear10 cells derived from a patient;

15

20

25

sample;

- (b) determining the concentration of Fc-receptor positive accessory cells in said first sample, and if required, supplementing said first sample with a second sample of mononuclear cells derived from said patient to achieve a concentration of 0.1-50% Fc-receptor positive accessory cells in said first sample;
- (c) contacting said first sample with OKT3 at or below 37°C to produce a T3CS;
 - (d) removing said T3CS from said first
- (e) determining the concentration of OKT3 in said T3CS, and if required, supplementing said T3CS with said OKT3 to achieve a concentration of at least 0.1 ng/ml;
- (f) providing a third sample of mononuclear
 cells derived from said patient;
 - (g) inactivating suppressor cells in said third sample; and
- (h) contacting said third sample with said 30 T3CS and said OKT3 for a period of time sufficient to activate said third sample in vitro to yield a population of immunoreactive cells.

- 11. The process of claim 10, wherein said Fc-receptor positive cells are monocytes.
- 12. The process of claim 10, wherein said concentration of Fc-receptor positive accessory cells is 1-30%.

5

- 13. The process of claim 12, wherein said concentration of Fc-receptor positive accessory cells is 5-15%.
- 14. The process of claim 13, wherein said
 10 concentration of Fc-receptor positive accessory cells is
 10%.
 - 15. A process of producing a population of immunoreactive cells comprising the steps of:
- (a) contacting a first sample of mononuclear : 15 cells derived from a patient with OKT3 at or below 37°C to produce a T3CS;
 - (b) removing said T3CS from said first sample;
- (c) determining the concentration of TNF α in 20 said T3CS, and if required, supplementing said T3CS with TNF α to achieve a concentration of at least 5 pg/ml;
 - (d) providing a second sample of mononuclear cells derived from said patient;
- (e) inactivating suppressor T-cells in said
 25 second sample; and
 - (f) contacting said second sample with said T3CS and for a period of time sufficient to activate said second sample *in vitro* to yield a population of immunoreactive cells.

- 45 -

- 16. The process of claim 15, wherein said TNF α concentration is in the range of 500-1000 ng/ml.
- 17. The process of claim 15, wherein step (a) is carried out for 6 hours-7 days.
- 5 18. The process of claim 17, wherein step (a) is carried out for 3 days.
 - 19. A process of producing a population of immunoreactive cells comprising the steps of:
- (a) contacting a first sample of mononuclear 10 cells derived from a patient with OKT3 at or below 37°C to produce a T3CS;
 - (b) removing said T3CS from said first sample;
- (c) determining the concentration of interferon-gamma (IFN-γ), interleukin-4 (IL-4), and interleukin-10 (IL-10) in said T3CS, and if required, adjusting the concentration of IFN-γ to at least 500 pg/ml, the concentration of IL-4 to less than or equal to 20 pg/ml, and the concentration of IL-10 to less than or equal to 20 pg/ml;
 - (d) providing a second sample of mononuclear cells derived from said patient; and
- (e) contacting said second sample with said
 T3CS and for a period of time sufficient to activate said
 second sample in vitro to yield a population of immunoreactive cells.

- 20. A process of producing a population of immunoreactive cells comprising the steps of:
- (a) contacting a first sample of mononuclear cells derived from a patient with OKT3 at or below 37°C to produce a T3CS;
- (b) removing said T3CS from said first sample;

10

20

- (c) determining the concentration of TNFα in said T3CS, and if required, supplementing said T3CS with TNFα to achieve a concentration of at least 5 pg/ml and determining the concentration of OKT3 in said T3CS, and if required, supplementing said T3CS with said OKT3 to achieve a concentration of at least 0.1 ng/ml;
- (d) providing a second sample of mononuclear
 15 cells derived from said patient;
 - (e) inactivating suppressor T-cells in said second sample; and
 - (f) contacting said second sample with said T3CS for a period of time sufficient to activate said mononuclear cells in vitro to yield a population of immunoreactive cells.
 - 21. The process of claim 6, further comprising after step (e):
- (f) measuring the number of CD25-positive cells in said second sample and discarding said second sample if said number is less than 10% of the total number of T cells in said second sample.
 - 22. The process of claim 6, further comprising after step (e):
- 30 (f) measuring the level of activation of said second sample and discarding said second sample if said level is less than twice the level of activation of a sample of unprocessed mononuclear cells.

- 23. The process of claim 2, further comprising after step (f):
- (g) measuring the number of CD25-positive cells in said second sample and discarding said second sample if said number is less than 10% of the total number of T cells in said second sample.

5

- 24. The process of claim 2, further comprising after step (f):
- (g) measuring the level of activation of said second sample and discarding said second sample if said level is less than twice the level of activation of a sample of unprocessed mononuclear cells.
- 25. The process of claim 4, wherein said T3CS comprises interleukin-1-alpha (IL-1α), interleukin-1-beta (IL-1β), interleukin-6 (IL-6), interleukin-8 (IL-8), TNFα, tumor necrosis factor-beta (TNFβ), IFNγ, granulocyte macrophage-colony stimulating factor (GM-CSF), monocyte/macrophage colony stimulating factor (M-CSF), and OKT3.
- 26. The process of claim 25, wherein said T3CS comprises:

		Amount ±	40%
	IL-1α	105	pg/ml
	IL-1β	1433	pg/ml
25	IL-6	808	pg/ml
	IL-8	213	ng/ml
	TNFα	570	pg/ml
	$\mathtt{TNF}\boldsymbol{\beta}$	171	pg/ml
	IFNγ	14350	pg/ml
30	GM-CSF	840	pg/ml
	M-CSF	1193	pg/ml
	ОКТ 3	12.7	ng/ml.

- 48 -

27. The process of claim 25, wherein said T3CS comprises:

		Amount ±	10%
	IL-1α	105	pg/ml
5	IL-1β	1433	pg/ml
	IL-6	808	pg/ml
	IL-8	213	ng/ml
	TNFα	570	pg/ml
	$\mathtt{TNF}\boldsymbol{\beta}$	171	pg/ml
10	IFNγ	14350	pg/ml
	GM-CSF	840	pg/ml
	M-CSF	1193	pg/ml
	OKT3	12.7	ng/ml.

- 28. The process of claim 26 or 27, wherein said
 15 T3CS comprises less than 20 pg/ml of each of interleukin2, interleukin-3, interleukin-4, interleukin-7,
 interleukin-10, interleukin-12, T cell growth factorbeta, and granulocyte-colony stimulating factor.
- 29. The process of claim 26 or 27, wherein said 20 T3CS comprises less than 100 units/ml of interleukin-2.
 - 30. The process of claim 6, further comprising contacting said second sample with an antigen.
- 31. A process of producing a population of antigen-specific polyclonal T cells comprising the steps of:
 - (a) contacting a first sample of mononuclear cells derived from a patient with OKT3 at or below 37°C to produce a T3CS;
- (b) removing said T3CS from said first
 30 sample;

- (c) determining the concentration of OKT3 in said T3CS, and if required, supplementing said T3CS with said OKT3 to achieve a concentration of at least 0.1 ng/ml;
- 5 (d) providing a second sample of mononuclear cells derived from said patient;
 - (e) contacting said second sample with said T3CS and an antigen for a period of time sufficient to activate said second sample in vitro to yield a population of antigen-specific polyclonal T cells.

10

- 32. The process of claim 31, wherein said OKT3 concentration is in the range of 0.1-1 ng/ml.
- 33. The process of claim 31, wherein said period of time is 1-30 days.
- 34. The process of claim 31, wherein said antigen is a tumor antigen.
 - 35. The process of claim 31, wherein said antigen is a bacterial antigen, viral antigen, or autoantigen.
- 36. The process of claim 1, further comprising 20 after step (e), contacting said second sample with IL-2.
 - 37. The process of claim 36, wherein said second sample is contacted with said IL-2 at 4°C.
- 38. The process of claim 36, wherein said IL-2 is present in an amount sufficient to bind to at least 25% of the IL-2 receptors on the surface of said immunoreactive cells.

- 39. The process of claim 38, wherein said IL-2 is present in an amount sufficient to bind to 100% of the IL-2 receptors on the surface of said immunoreactive cells.
- 5 40. The process of claim 1, wherein step (e) is carried out below 37°C.
 - 41. A population of immunoreactive cells produced by the process of claim 1, wherein at least 75% of said cells are T lymphocytes.
- 42. An immunoreactive mononuclear cell produced by the process of claim 1.
 - 43. The cell of claim 42, wherein said mononuclear cell is a T cell.
- 44. The cell of claim 43, wherein said mononuclear cell is a Th1-type T cell.
 - 45. The cell of claim 43, wherein said cell expresses at least 10% more cell-surface CD25 than an unprocessed T cell.
- 46. The cell of claim 43, wherein said cell
 20 proliferates at a rate that is at least twice that of an
 unprocessed T cell when contacted with an immune
 stimulant.

- 51 -

- 47. A process for producing a mixture of autologous nonspecific lymphocyte activators comprising the steps of:
- (a) providing a sample of mononuclear5 cells derived from a patient; and
 - (b) contacting said sample with a CD3-binding compound at or below 37°C, wherein said CD3-binding compound is present at a concentration of at least 0.1 ng/ml.
- 10 48. A process for producing a mixture of autologous nonspecific lymphocyte activators comprising the steps of:
 - (a) providing a sample of mononuclear cells derived from a patient; and
- 15 (b) contacting said sample with a T cell receptor-binding compound at or below 37°C, wherein said T cell receptor-binding compound is present at a concentration of at least 0.1 ng/ml.
- 49. A mixture of autologous nonspecific
 20 lymphocyte activators produced by the process of claim
 47.
 - 50. The mixture of claim 49, wherein said CD3-binding compound is OKT3.
- 51. The mixture of claim 50, wherein the concentration of said OKT3 is at least 0.1 ng/ml.
 - 52. The mixture of claim 51, wherein said mixture comprises at least 5 pg/ml $TNF\alpha$.

- 52 -

53. The mixture of claim 51, wherein said mixture comprises:

		<u>Amount ±</u>	10%
	IL-1α	105	pg/ml
5	IL-1β	1433	pg/ml
	IL-6	808	pg/ml
	IL-8	213	ng/ml
	$\mathtt{TNF}\alpha$	570	pg/ml
	$\mathtt{TNF}\boldsymbol{\beta}$	171	pg/ml
10	IFNγ	14350	pg/ml
	GM-CSF	840	pg/ml
	M-CSF	1193	pg/ml
	OKT3	12.7	ng/ml.

54. The mixture of claim 51, wherein said mixture 15 comprises:

		Amount ±	40%
	IL-1α	105	pg/ml
	IL-1β	1433	pg/ml
	IL-6	808	pg/ml
20	IL-8	213	ng/ml
	TNFα	570	pg/ml
	$\mathtt{TNF}\boldsymbol{\beta}$	171	pg/ml
	IFNγ	14350	pg/ml
	GM-CSF	840	pg/ml
25	M-CSF	1193	pg/ml
	ОКТ3	12.7	ng/ml.

- 55. The process of claim 47, wherein step (b) is carried out for 1-7 days.
- 56. The process of claim 55, wherein step (b) is 30 carried out for 3 days.

- 57. The process of claim 47, wherein said CD3-binding compound is OKT3.
- 58. The process of claim 47, further comprising contacting said sample with an antigen.
- 5 59. The process of claim 47, wherein step (b) is carried out at 29-34°C.
 - 60. A pharmaceutical composition comprising an immunoreactive cell produced by the process of claim 1 in a pharmaceutically acceptable carrier or diluent.
- 10 61. A method of treating a malignant tumor in a patient comprising administering to said patient an immunoreactive mononuclear cell produced by the process of claim 1.
- 62. The method of claim 61, wherein a suppressor cell-inhibiting compound is concurrently administered to said patient.
 - 63. The method of claim 62, wherein said suppressor cell-inhibiting compound is cimetidine, indomethacin, or both cimetidine and indomethacin.
- 20 64. The method of claim 61, wherein said tumor is selected from a group consisting of renal cell carcinoma, breast carcinoma, prostate carcinoma, colo-rectal carcinoma, pancreatic carcinoma, ovarian carcinoma, melanoma and non-small cell carcinoma of the lung.
- 25 65. The method of claim 61, wherein said tumor is a leukemia or lymphoma.


- 54 -

66. The method of claim 64, wherein said tumor is a renal cell carcinoma.


- 67. A process of producing a population of immunoreactive cells comprising the steps of:
- 5 (a) contacting a first sample of mononuclear cells derived from a patient infected with a pathogenic virus with OKT3 at or below 37°C to produce a T3CS;
 - (b) removing said T3CS from said first
 sample;
- (c) determining the concentration of TNFα in said T3CS, and if required, supplementing said T3CS with TNFα to achieve concentration of at least 5 pg/ml or determining the concentration of OKT3 in said T3CS, and if required, removing said T3CS to achieve concentration less than 0.1 ng/ml;
 - (d) providing a second sample of mononuclear cells derived from said patient;
 - (e) inhibiting suppressor T-cells from said second sample; and
- T3CS for a period of time sufficient to activate said second sample of mononuclear cells in vitro to yield a population of immunoreactive cells.
- 68. The process of claim 67, further comprising 25 after step (f):
 - (g) measuring the level of activation of said second sample and discarding said second sample if said level is less than twice the level of activation of a sample of unprocessed mononuclear cells.


- 55 -

- 69. The process of claim 67, wherein said virus is hepatitis B virus, hepatitis C virus, herpes simplex virus, varicella-zoster virus, cytomegalovirus, papilloma virus, human immunodeficiency virus-1, human immunodeficiency virus-2, or Epstein Bar virus.
 - 70. An immunoreactive mononuclear cell produced by the process of claim 67.
- 71. A pharmaceutical composition comprising a mononuclear cell produced by the process of claim 67 in a pharmaceutically acceptable carrier or diluent.
 - 72. A method of treating a patient infected with a pathogenic virus comprising administering to said patient the cell of claim 70.
- 73. The method of claim 72, wherein said pathogenic virus is hepatitis B virus.
 - 74. The process of claim 67, wherein said patient is infected with hepatitis B virus.
 - 75. The process of claim 31, wherein said antigen is a hepatitis B-derived antigen.

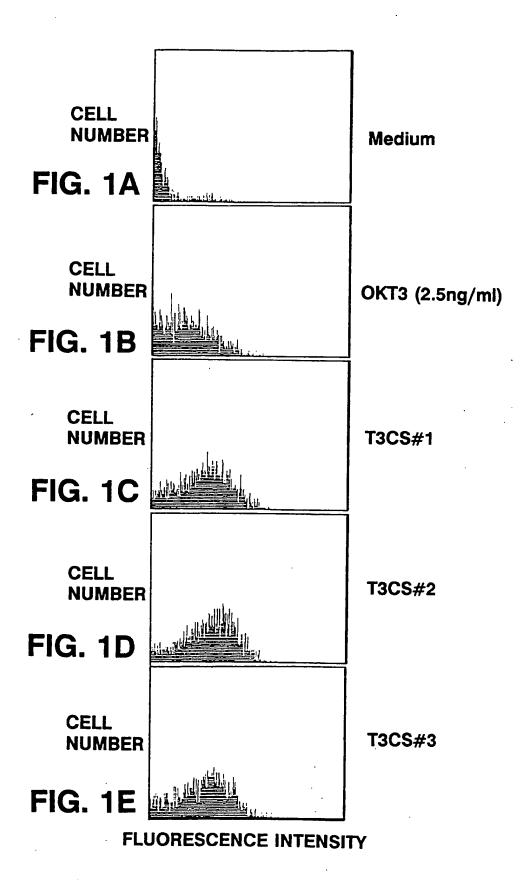



SUBSTITUTE SHEET (RULE 26)

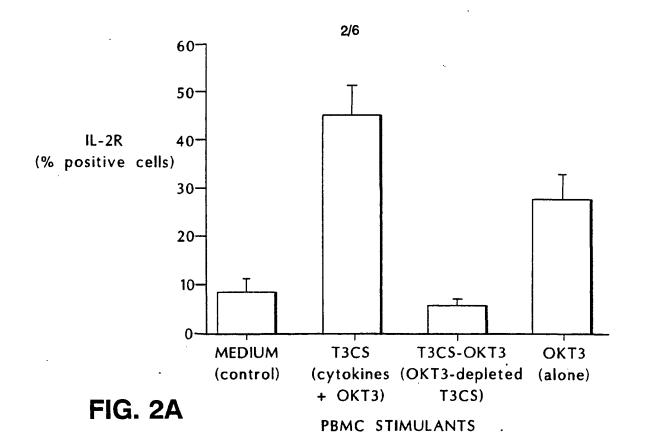


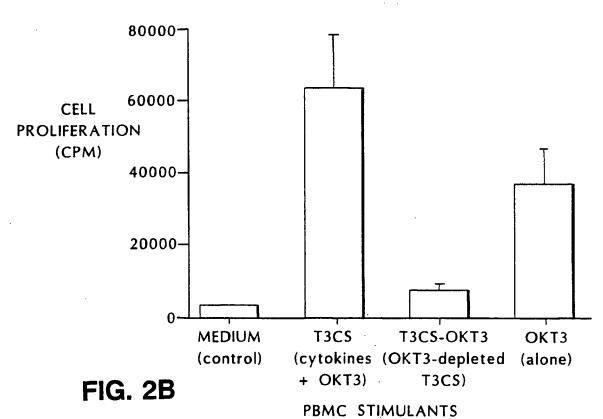


SUBSTITUTE SHEET (RULE 26)

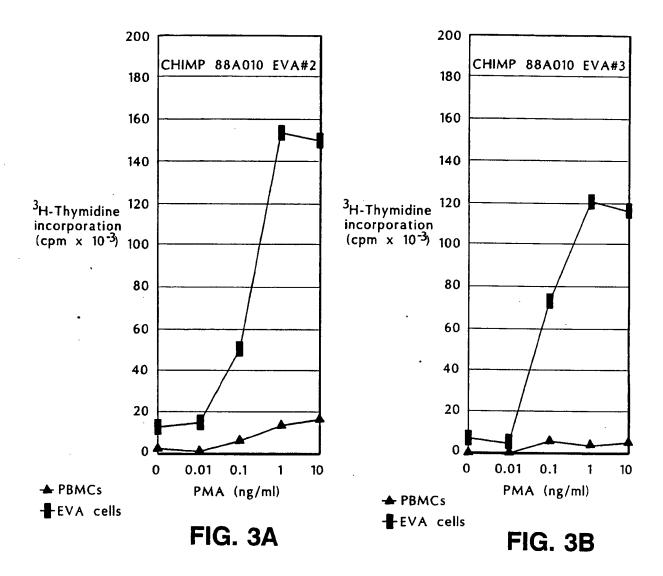

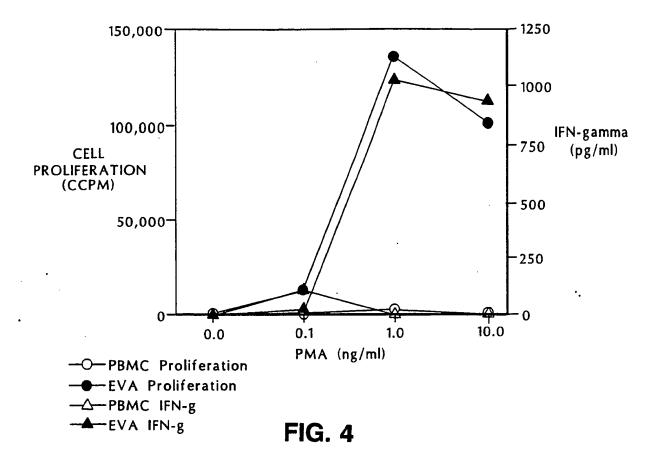

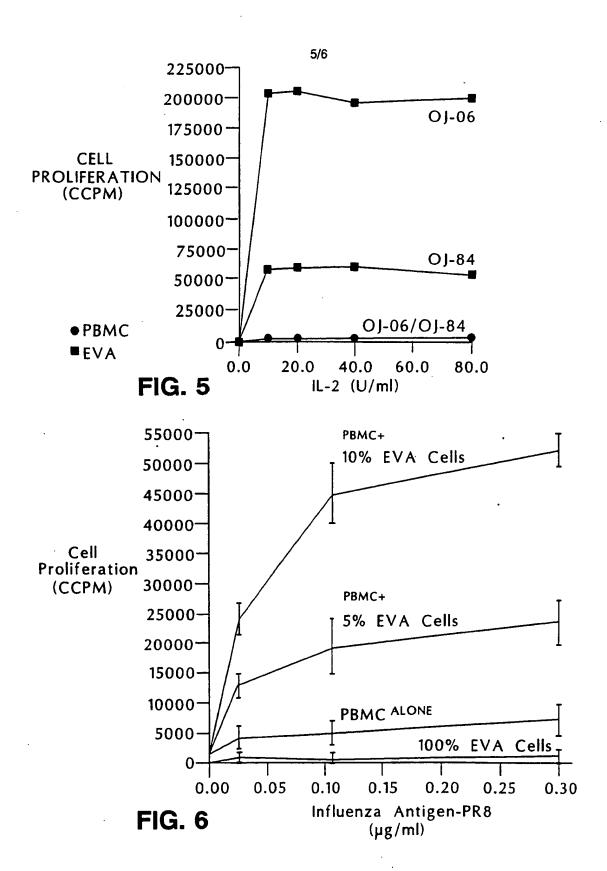

6/6

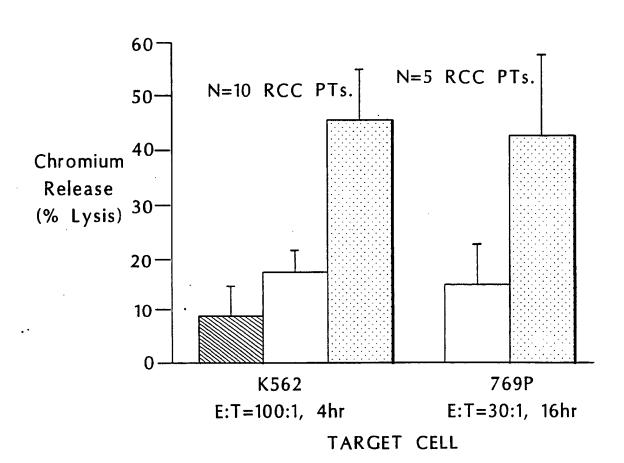

PBMC

Medium Control

K562= Leukemia line 769P= RCC cell line

EVA Cells

FIG. 7

INTERNATIONAL SEARCH REPORT

Incomational application No. PCT/US95/00193

A. CLASSIFICATION OF SUBJECT MATTER			
IPC(6) :C12N 5/08; A61K 121/00			
	:435/240.2, 240.21; 424/93.1, 93.21		
	to International Patent Classification (IPC) or to both	national classification and IPC	
	LDS SEARCHED		
Minimum d	locumentation searched (classification system followed	d by classification symbols)	
U.S. :	435/240.2, 240.21; 424/93.1, 93.21		
Documenta	tion searched other than minimum documentation to th	e extent that such documents are included	l in the fields searched
Electronic o	data base consulted during the international search (n	ame of data base and, where practicable	search terms used)
i .	i, APS, EMBASE, WPI, LIFESCI, MEDLINE, BIO		,
C. DOC	CUMENTS CONSIDERED TO BE RELEVANT		
Category*	Citation of document, with indication, where a	ppropriate, of the relevant passages	Relevant to claim No.
Y	TRANSPLANTATION PROCEEDINGS, VOL. 24, NO.6, ISSUED 01 DECEMBER 1992, P.T. LAVIN ET AL., "AUTOLYMPHOCYTE THERAPY FOR METASTATIC RENAL CELL CARCINOMA: INITIAL CLINICAL RESULTS FROM 335 PATIENTS TREATED IN A MULTISITE CLINICAL PRACTICE", PAGES 3059-3064, ENTIRE DOCUMENT.		1-75
Y	US, A, 5,192,537 (OSBAND) 09 MARCH 1993, COLUMS 1- 10, SEE ENTIRE DOCUMENT.		1-75 ·
Y	US, A, 4,716,111 (OSBAND ET A		1-75
Further documents are listed in the continuation of Box C. See patent family annex.			
l	ocial categories of cited documents:	"I" inter document published after the inte date and not in conflict with the applier	ation but cited to understand the
	cument defining the general state of the art which is not considered be of particular relevance	principle or theory underlying the inve	ention
'E' cas	rlier document published on or after the international filing date	"X" document of particular relevance; the considered novel or cannot be consider	e claimed invention cannot be red to involve an inventive step
	cument which may throw doubts on priority claim(s) or which is and to establish the publication date of another citation or other	when the document is taken alone	,
	ecial remon (as specified)	"Y" document of particular relevance; the considered to involve an inventive	
	cument referring to an oral disclosure, use, exhibition or other	combined with one or more other such being obvious to a person skilled in th	documents, such combination
"P" document published prior to the international filing date but later than "&" document member of the same patent family the priority date claimed			
Date of the actual completion of the international search 13 FEBRUARY 1995 Date of mailing of the international search 0 9 MAR 1995			
Commissioner of Patents and Trademarks		Authorized officer Dilly No	There I'm
Box PCT Washington, D.C. 20231 CHRISTOPHER EISENSCHENK			
1	in (703) 305-3230	Telephone No. (703) 308 0106	